

A close-up, macro photograph of a woven brush mat. The mat is made of dark brown, textured fibers, likely horsehair or a similar natural material, woven in a tight, regular pattern. The lighting is dramatic, highlighting the texture and creating a sense of depth. The background is blurred, focusing attention on the intricate weave.

M.P.F.

BRUSH COMPANY

2011
2012



- Introduction 04
- Brushes and Brush Kits 06
- Accessories 14
- Palettes 24
- Evolution Series 32
- Wax and Die Management 40
- Rotary Instruments 52
- Techlight 66
- Pilot Labs 72

table of contents



Ladies & Gentlemen,

High quality functional products, reasonable pricing, and attention to detail are the key ingredients for a successful product range. We at **M.P.F. Brush Company** put a lot of energy into these basic factors so that you receive a high quality product at an affordable price.

Our primary goal is to continue to improve our customer's performance by offering a unique range of laboratory hand tools, brushes, palettes, accessories, rotary instruments, wax and die management systems and dental lab LED bench lighting.

We are very proud of our latest innovation the **M.P.F. Evolution series**, we have combined materials such as aluminum, carbon fiber & resin to create this unique range of dental ceramic hand tools.

We feel all technicians worldwide will appreciate, the attention to detail of this exclusive series and its form and function which makes this new range unique.

Our collaboration with opinion leaders from around the world has allowed us to "think outside the box" when designing unique dental ceramic hand tools.

The rapid growth of the **M.P.F. product line** worldwide is something we highly appreciate and we would like to take this opportunity to thank all our customers from around the world.

President
George Papadopoulos.
CDT

**Committed to
excellence, quality
& innovation**

10/10

**Feedback from
our users means
we're getting it right**

snapping kolinsky

Brushes & Brush Kits



Brushes & Brush Kits

The unique spring activated 'flick' action of the **M.P.F. Ceramic Brush**, reshapes and contours the brush hair with one simple movement. This 'flick' action provides the perfect tip ready for application, time after time. The versatile, replaceable, 'spring activated head', allows the dental ceramist fantastic manoeuvrability. It's patented 'spring activated head' makes it possible to use this brush in ways no other conventional brush has ever been used before. The 'flick' action also allows the dental ceramist greater ease when cleaning the brush tips. The 'spring activated head' makes cleaning effortless and helps to maintain the integrity of the brush. Simply apply gentle pressure to the head of the brush using your fore-finger pushing it towards you to release the spring, it couldn't be easier. This action expels the excess water and paste from the brush without wiping or handling of the brush hair in any way, while at the same time, creating the finest point you have ever seen on a brush tip...every time!

The M. P. F. Spring Activated Ceramic Brush handles are made from lightweight aluminum. They are anodized in 6 contemporary colors, red, blue, black, gold, chrome and purple.

We recently added a **NEW SERIES** of **RESIN** handles to our collection. These resin handles are light weight, strikingly designed and available in 5 different colors, red, blue, green, black and orange.

Our unique 'Replaceable Brush Tip System' has the advantage over all other brush systems. Our handles are made to fit all size brush tips. This is due to the unique design of our screw on system, which attaches itself to the brush tip of your choice. Our series of 3 different sized replaceable brush tips, 4, 6 and 8, are available in either the finest 100% male Siberian Kolinsky hair, or fine synthetic hair.

A very popular item in our catalogue is the **3 in 1 Brush Kit**. Our kit provides the ceramist with **a complete set of three brush tip sizes, and one versatile handle that fits all.**

All our brush tips are made in Germany by the finest brush makers in the world!

The M.P.F. Brush Company has a brush for every technician, and a complete range of wooden handled brushes for a more economical solution!

Spring activated ceramic brush
100%
moisture control
with one simple flick

SPRING ACTIVATED CERAMIC BRUSHES

Aluminum Handle



Resin Handle



Colors Available

Aluminum Anodized Handles


Available in 6 brilliant colors:



Colors Available

Resin Handles

Available in 5 brilliant colors:

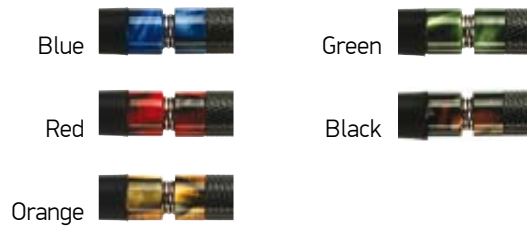

NON SPRING ACTIVATED CERAMIC BRUSHES

Aluminum Handle




Resin Handle



Colors Available

Aluminum Anodized Handles

Available in 6 brilliant colors:

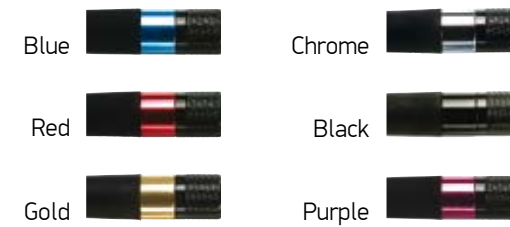

Colors Available

Resin Handles

Available in 5 brilliant colors:



BRUSH KITS

Aluminum Handle



Spring Activated
Item Order No: 004
Available in all 6 colors



Non Spring Activated
Item Order No: 012
Available in all 6 colors

Resin Handle




Spring Activated
Item Order No: 008
Available in all 5 colors



Non Spring Activated
Item Order No: 016
Available in all 5 colors

REPLACEMENT BRUSH TIPS

100% Kolinsky Hair

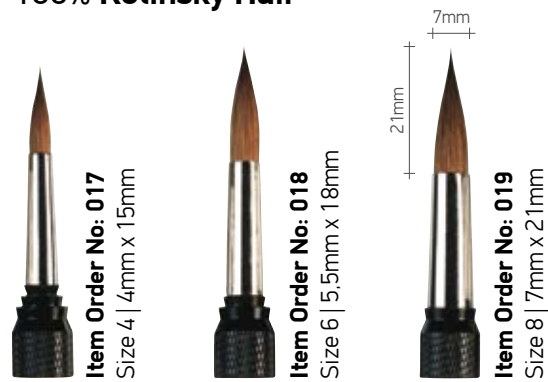

100% Synthetic Hair



Stain Brush



Item Order No: 023

The most important aspect of a good Stain brush is the quality of the hair and the sharpness of the tip. Our stain brush meets both these requirements. Our Stain brush is manufactured with the highest quality Kolinsky hair, creating a perfectly designed sharp brush tip.



Opaque Brush



Item Order No: 024

Our Opaque brush has a unique 'cats tongue' brush tip design, allowing large amounts of opaque to be applied. Its very sharp point allows the ceramist flexibility during detailed applications. Manufactured with the highest quality Kolinsky hair.



Contour Brush



Item Order No: 025

Our Contour brush has a unique fan shaped head, allowing the ceramist to contour and shape semi dry porcelain whilst the ceramic build up is taking place. The hair is a mix of kolinsky and very fine synthetic hair. Such a hair combination leaves no striations on the ceramic build up. The fan shaped head allows access to the finest embrasures without damaging the detailed build up.



WOODEN HANDLED CERAMIC BRUSHES

Occlusal Brush



Item Order No: 026


The Occlusal brush is a very small sized contour brush. During detailed Occlusal design we have found that this size brush head is wonderful for contouring and shaping intricate detailed occlusal surfaces. Like the contour brush it is made from a combination of Kolinsky and synthetic hair.



Classic Brushes



Item Order No: 027 Classic Brush 4

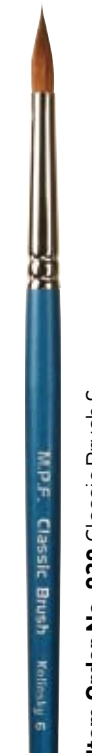

Item Order No: 028 Classic Brush 6



Item Order No: 029 Classic Brush 8

Our M.P.F. Classic brushes are produced using the highest quality Kolinsky hair. Inside these brush heads are several metallic beads, which act as a condenser when shaken. This action sends the moisture trapped within the brush head, directly to the brush tip, giving a very sharp point to work with. The handles are made from wood, varnished in a matt blue color. The Classic brushes come in the traditional three sizes 4, 6 & 8.

precise
functional
unique

Accessories



Accessories

Manufacturing dental ceramics throughout the ages has relatively been unchanged as far as the manufacturing process is concerned. What has changed immensely are the tools we use to manufacture great ceramic restorations.

At the **M.P.F. Brush Company** we strive to create dental ceramic accessories that not only look great, but their functionality is a great help to the ceramist's daily grind at the bench.

Our collaboration with our pilot lab EU and pilot lab USA, we have established an innovative design team. This team effort has created some of the finest dental ceramist accessories the market has seen. Utilizing materials such as aluminum, resin, carbon fiber & wood we have had the ability to create attractive, functional products, and all at an affordable price.

Dental laboratory hand tools designed and made by technicians, for the technician!!!

Making your work
as easy as

123

made by technicians
for technicians

MASTER BRUSH KIT



1. Master Opaque Brush
2. Master Occlusal Brush
3. Master Spring Activated Brush. Size 8
4. Master Stain Brush
5. Master Contour Brush

My collaboration with the M.P.F. Brush Company this past year has given me the ability to create this new personalized product line. Utilizing light weight Aluminum and the purest Mahogany wood, we have created some innovative dental ceramic hand tools. The Master brush kit gives you a range of brush shapes that are unique in design and excellent in hair quality. The Master Ceramic brush is made with the highest quality Siberian Kolinsky hair, my personal input into the design of the brush tip and what I feel the dental ceramist needs most. A fine sharp Kolinsky brush tip, with a strong belly of hair in the mid section for extra leverage, so that you can pick up any amount of ceramic material you need, but also the finest point for detailed anatomical applications eg. Occlusal Area. The fine sharp point of the Master stain brush allows you to create the finest of details, the Opaque brush with its cat tongue design gives you a nice sharp point, but also the ability to apply larger portions of opaque when needed. The Contour brush and Occlusal brushes are made with the finest synthetic hair, this hair is soft yet firm, allowing the technician the ability to shape and contour detailed anatomy into the ceramic build up.

- Complete Master Brush Kit is available either,**
 With size 4 Kolinsky Brush Tip **Item Order No: 065A**
 With size 6 Kolinsky Brush Tip **Item Order No: 065B**
 With size 8 Kolinsky Brush Tip **Item Order No: 065C**
 Replacement Master Contour Brush Tip **Item Order No: 066**
 Replacement Master Stain Brush Tip **Item Order No: 067**
 Replacement Master Opaque Brush Tip **Item Order No: 068**
 Replacement Master Occlusal Brush Tip **Item Order No: 069**

- All components of the Master Brush kit can be purchased individually**
 Complete Master Brush + handle size 4 brush tip **Item No: 072**
 Complete Master Brush + handle size 6 brush tip **Item No: 073**
 Complete Master Brush + handle size 8 brush tip **Item No: 074**
 Complete Master Brush + handle Contour brush **Item No: 078**
 Complete Master Brush + handle Stain brush **Item No: 079**
 Complete Master Brush + handle Opaque brush **Item No: 080**
 Complete Master Brush + handle Occlusal brush **Item No: 081**

- All components of the Master Carver kit can be purchased individually**
 Complete Master Carver handle flat blade **Item No: 082**
 Complete Master Carver handle Spear **Item No: 083**
 Complete Master Carver handle Contra Angle **Item No: 084**
 Complete Master Carver handle Dagger **Item No: 085**

MASTER CARVER



Complete Master Carver
Item Order No: 048



1. "Contra Angle" Replacement Blade. **Item Order No: 054**
2. "Spear" Replacement Blade. **Item Order No: 053**
3. "Dagger" Replacement Blade. **Item Order No: 055**
4. "Flat" Replacement Blade. **Item Order No: 052**
5. Stainless Steel Master "Klip" Tweezers. **Item Order No: 056**

The Master Carver again has been made utilizing the same materials Aluminum and Mahogany wood, but also uses surgical stainless steel for the replaceable blades.

I have designed four shapes which I feel cover all the needs for the dental ceramist .



The Master Klip set of tweezers has been designed to have enough grip on the finest veneers or all ceramic crowns, without causing any damage to your delicate pieces of work.

The combination of Aluminum and Mahogany wood has also been used for the stands in both product lines. Attention has been given to the delicate size of the complete range, so that the ceramists table is not over crowded with bulky instrument stands.

A must for me was that all instruments had replaceable brush tips & replaceable blades, after our everyday grind, instruments sustain some wear & tear, so replacement parts are a must.

I hope that you find the Master Series fresh & innovative.

Sincerely,
Nondas Vlachopoulos. CDT

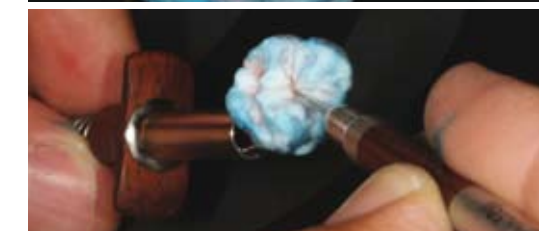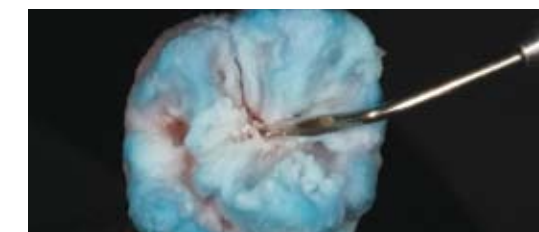

CERAMIC BRUSH HYDRATOR



Item Order No: 035

This new innovative idea is a must for all dental ceramists. Kolinsky hair brushes tend to be difficult to moisten and prepare for build up, so we've designed a water reservoir that suspends the brush without touching the sides or bottom of the vessel. After five minutes of soaking, the water slowly permeates its way into the hair bundle, giving the perfect moisture level required for application.

NB: The hydrator is "not" to be used as an over night stand to suspend the brush. If the brush is left for a long period of time in water, problems to the hair may occur. The water in the hydrator should be changed on a regular basis. **NB:** Brush not included

CERAMIC TABLE MAT



Item Order No: 042
35cm x 55cm

The M.P.F. Ceramic work bench mat is a 1.5mm thick double layered plastic sheet which is 55cm in width and 35cm in height. This wonderfully designed item displays the dental art of dental ceramist Nondas Vlachopoulos.

This assortment of hand made ceramic teeth is an exceptional piece of dental art to say the least! A wonderful accessory for the dental ceramist!

HONEYCOMB CROWN & BRIDGE STAND



Item Order No: 059
Honey Comb Ceramic CR & BR Stand.
With 10 Stainless Steel Pins and 10 Zirconia Pins.



Item Order No: 036
Honey Comb Ceramic CR & BR Stand.
With 20 Stainless Steel Pins.



Item Order No: 051
Honey Comb Ceramic CR & BR Stand.
With 20 Zirconia Pins.

This Honeycomb CR & BR ceramic firing tray is a ceramic based material which can withstand temperatures of up to 1500 degrees Celsius.

NYLON CERAMIC PALETTE WORKSHEET



Item Order No: 032

Replacement Ceramic Palette Nylon work sheet. 2 sheets per pack.

BRUSH POOL



M.P.F. Brush Pool Item Order No: 043
Brush not included



M.P.F. Replacement Blue Plastic Water Bowl
Item Order No: 044



M.P.F. Brush Pool Replacement Sponge
Item Order No: 045

The M.P.F. Brush Pool was designed to complement the M.P.F. spring activated brush system. This Unique design incorporates a water bowl, sponge and protective backboard all in the one unit. When dipping the brush into the bowl to clean, we must then flick the excess water from the brush, by performing this action, excess water is driven out of the brush head, which often creates a messy work space, now with the M.P.F. Brush Pool, the excess water is contained against the backboard, preventing the excess water from creating a messy work bench.

The M.P.F. Brush Pool is made from very durable plastic and can be washed or steamed without any problem.

The M.P.F. Brush Pool has been designed to suit both left and right handed technicians!



AGATE SPATULA



Item order No: 030

Our agate spatula is made from raw Brazilian Agate Stone, shaped in such a way to assist the ceramist when mixing ceramic powder.

The tail end of the spatula is also shaped at an angle, so small amounts of ceramic powder can be taken from bottles with ease.

NB: The agate spatula being a natural material may vary in color.



DIAMOND PASTE



Item Order No: 058

The M.P.F. diamond paste should be used after the ceramic or composite surface has been treated with our fine rubber polishers. We then collect small amounts of paste onto a felt mop, and then at 8000 -12000 rpm applying slight pressure, a high shine luster will appear.

Each Jar Contains 11gr of Diamond Paste.



SPIDER KLIP



Item Order No: 071

No longer do you need to use 3-4 tweezers to stain your single anteriors! By simply loading all your crowns or veneers into the Spider Klip tweezers, you have a perfect arched visual and your staining process can be precise, due to the fact that you have all the crowns lined up in front of you. With the six spring activated tweezers all in one simple handle, the ceramist has a much better visual of the crowns to be stained.

To stain 6 anterior crowns or veneers, with the conventional method, you must use at least 2 tweezers and constantly pick up one crown, position it and stain it, then leave it down and pick up another set of tweezers and repeat the process. This is not only time consuming, but also create an inaccurate color staining process.



1.



2.



3.



4.

1. Six stainless steel spring activated tweezers aligned into the shape of the dental arch.
2. By activating the spring, tweezers open, to give you enough room to position your crown or veneer. The rubber grip also acts as a dampener, to prevent any rotation of your veneer or crown while the staining procedure takes place.
3. Position crown in tweezers: 6 crowns can be positioned in less than 30 seconds!
4. Crowns positioned ready to proceed with staining.

This latest innovation from the M.P.F. Brush Company, has been designed to allow the dental ceramist a much faster and more accurate way to stain up to 6 crowns or veneers at the one time, thus making the staining process much more precise and color correct.

With the M.P.F. Spider Klip, it can take you less than 4 minutes to load 6 crowns or veneers into the hand piece and complete your staining procedure. When staining on the M.P.F. Spider Klip, the process is far more accurate, due to the fact the crowns are lined up next to each other, replicating the arch form of the dentition.

The Spider Klip is made of aluminum and the tweezers and spring are stainless steel.

A soft rubber grip has been positioned on the posterior of the tweezers, acting as a dampener and preventing the crown or veneer from moving around while applying your stain material.

The front tweezers are shaped in such a way that there are two points of contact against the anterior surface of the crown or veneer. This combination of rubber grip and 3 point contact make for a very stable positioning of the item to be stained.

The Spider Klip also comes with 6 replacement rubber grips.

If you wish to work without the rubber grips positioned on the tweezers, they can be removed.

humid
consistent
reliable

Palettes



Palettes

Palettes are a tool utilized by the dental ceramist every day, and are a necessary accessory in the ceramic laboratory.

As a tool, they must firstly be effective, simple to use and assist the technician to be productive. Secondly, they must be aesthetically pleasing; as it is likely the technician will be looking at it for a good part of the working day.

With its combination of natural stone and durable plastic base, the M. P. F. Ceramic Palette is very easy to work with. The simple 'nylon palette work sheet' allows the water from the base to permeate through the stone and naturally keep the ceramic mix constantly wet 24/7.

We have recently added a NEW edition to our ceramic palette range, the M.P.F. Glass Plate Ceramic Palette. This new palette is designed to cater to the ceramist that prefers to work on a glass plate, but at the same time have the humid effect utilizing the synthetic leather strips which permeate moisture from the base of the palette.

The Paste Opaque and Stain Palettes have a very unique design style which makes access to the ceramic medium, a breeze. The combination of 'durable plastic base', 'sintered ceramic reservoirs' and a 'tightly sealing lid' make our palettes an essential part of any busy laboratory.

Check out our latest ceramic palette innovation the M.P.F. Evolution Palette, you will be amazed!

Choose between
natural stone
or glass for humid
consistency

24/7

GLASS PLATE CERAMIC PALETTE



Item Order No: 087

The M.P.F. glass plate ceramic palette is one of the latest innovations by the M.P.F. Brush Company.

The M.P.F. glass plate ceramic palette consists of a glass plate and plastic base. It is supplied with one agate spatula for mixing ceramic, four synthetic leather strips and a white plastic bottle for distilled water. We have created this ceramic palette for the ceramists that prefer to work on a glass surface. **By utilising the synthetic leather strips** we have been able to create a glass plate ceramic palette that can keep your ceramic mix continuously wet. The synthetic leather strips are positioned across the glass plate in 4mm location grooves which lock the strip in position and **by placing your ceramic material up against these synthetic leather strips, the water from the base permeates** up through the strips and into the ceramic stack, creating a very homogenous ceramic mix. To obtain an even greater consistency special modeling liquids can be applied directly to the ceramic stack itself.

1. Replacement Leather Strips pack of 4.

Item Order No: 088

2. Replacement Water Bottle. **Item Order No: 089**

3. Replacement Glass Plate only. **Item Order No: 090**



1. Take the synthetic leather strips and soak them in water.

2. On the glass plate there are two 4mm grooves that have been scored into the glass, so that the synthetic leather strips can be positioned.


3. You will notice that the synthetic leather strips are slightly longer on both sides of the glass plate, this is so they can be rolled down into the plastic base to permeate the water from the base up into the ceramic stack.

4. Using the agate spatula provided position the synthetic leather strips down into the plastic base.

5. On the top right hand corner of the glass plate there is a blue metallic ringlet, we have created an entrance through the glass plate via this blue ringlet, which makes it easy to fill the base with distilled water, utilising the 250ml water bottle provided.

6. The M.P.F. Glass Plate Ceramic palette system, creates a very homogenous ceramic stack, if the ceramist desires a more homogenous feel to their ceramic, special modeling liquids can also be applied.

CERAMIC PALETTE

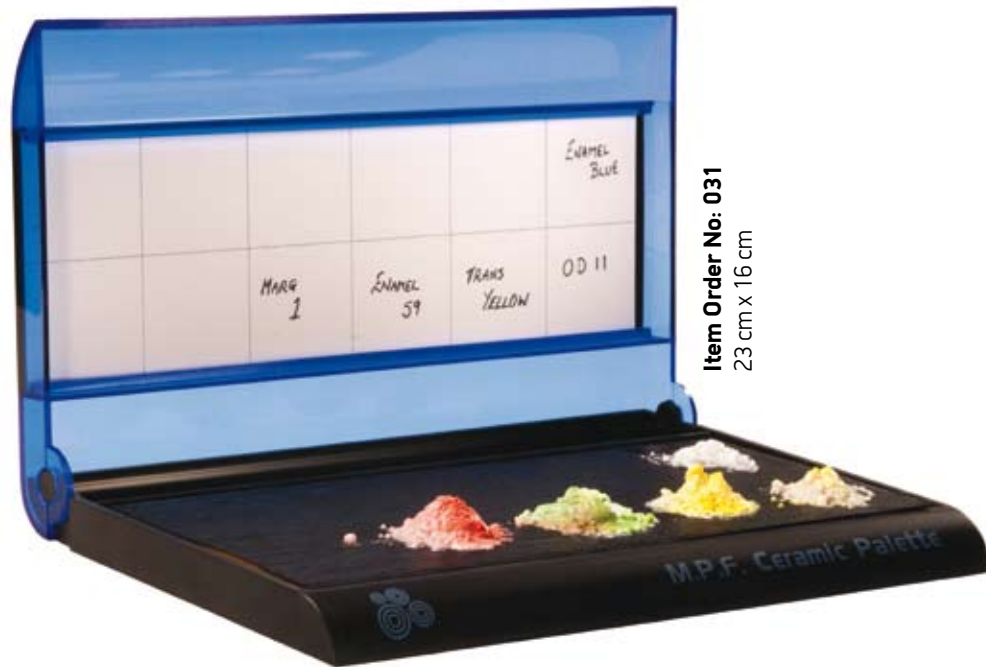

Item Order No: 031
23 cm x 16 cm

The unique marbel slab allows controlled amounts of water to permeate through the nylon worksheet.

M.P.F. Brush Company introduces a new advancement for the discriminating ceramist. The new M.P.F. Ceramic Palette is the latest innovation from the M.P.F. Brush Company. The palette comes in a black case with an Italian marble stone, using capillary action to wet ceramic powders.

This unique marble allows controlled amounts of water to permeate and wet the nylon worksheet, to allow uniform wetting of the ceramic. This will not allow over wetting of the ceramic material but provides the ceramist with a more consistent and workable ceramic mix. The nylon work sheet that covers the marble is kinder to Kolinsky brushes, like the M.P.F. Spring Activated Kolinsky Brush and helps to keep a finer point for longer periods of time. The palette comes with 2 nylon work sheets that are durable and easy to clean.

We also supply a plastic Shade palette work sheet, located on the internal side of the clear plastic lid. This sheet enables the ceramist to record the specific shade used on the palette at the time of each build up, using a permanent marker pen, which can be simply wiped clean with white spirit in order to be reused.

Nylon palette worksheets are also sold seperately. Item order No: 032.



1. Our Italian marble slab has natural pores of approximately 500 microns.
2. Our plastic shade palette work sheet is divided into quadrant. With a marker pen we can mark down a memory of what ceramic colours we are using for each build up.
3. Detail of ceramic mix.

STAIN & PASTE OPAQUE PALETTE



Item Order No: 033
Paste Opaque Palette 18 cm x 16 cm



Item Order No: 034
Stain Palette 18 cm x 16 cm

"The unique angled step type design, allows easy access to materials..."



Ceramic mixing cup
Due to the shape of the ceramic cup, all the material flows to the back.



1. The Ceramic cup can hold a considerable quantity of stain material.
2. The ceramic mixing pad allows us to experiment with several color combinations.
3. Individual cups allow us to control our color opaques and not waste material.

Our Stain & Paste Opaque palettes offer a great solution for storing stain and paste opaque materials. Both palettes have the designated positions for each color screen printed onto the plastic base so that the colors are kept in order.

The unique angled step type design, allows easy access to materials. The palette is constructed from durable plastic and the white reservoirs are made from glazed ceramic. The reservoirs have a very unique design to them, they are elongated towards the rear so that the material will flow to the back of the reservoir, this prevents any spillage.

The lid is made from a clear plastic and has been designed to create a tight seal around the black base, preventing too much oxygen entering the palette, which can cause the materials to dry out.

NB: When your palette is not in use always keep the lid closed, preventing the material from drying out and dust contamination.

f
form.
f
function
f
future

Evolution Series



Item Order No: 2025 M.P.F. Evolution Master Ceramic Palette replacement sponge

Item Order No: 2026 M.P.F. Evolution Master Ceramic Palette glass tablet

Item Order No: 2027 M.P.F. Evolution Master Ceramic synthetic absorbent strips pack of 4

Item Order No: 2000 M.P.F. Complete Evolution Master Ceramic Palette 31 x 18 cm

Evolution Series

At the M.P.F. Brush Company esthetics, form and function are the key elements when designing a new product.

The New M.P.F. Master Evolution Ceramic Palette has been developed taking all these elements into consideration, its ergonomic design incorporates the sponge for cleaning your brush, the slide drawer with its internal glass stain palette, and the unique synthetic leather strip which permeates the moisture from the reservoir, up into your ceramic mix, all this in one unit, making it neat and efficient. Utilizing the smooth glass plate tablet, we prevent our very delicate Kolinsky brushes from being destroyed.

For years the M.P.F. Brush Company has been developing precision instruments, designed by technicians for technicians, this High End M.P.F. Evolution series is no different, having the inspiration of Peter Pizzi M.D.T. we have combined our innovative ideas, and this stylish also functional series of products has evolved.

Our carbon fibre evolution series

100%

more efficient
more effective

EVOLUTION MASTER BRUSH KIT



The M.P.F. Master Evolution Brush assortment kit.

We have created this unique stylish brush assortment kit, combining a slick Aluminum handle with a mat carbon fiber grip for a comfortable feel and grip.

The Evolution brush assortment kit comes with 1 master brush with the brush tip of your choice, 4, 6 or 8, 1 Stain brush, 1 Contour brush for shaping semi dry ceramic in the buildup stage, and 1 flat glaze brush. All brush tips are replaceable and are produced with the highest quality Kolinsky hair.

Item Order No: 2001 M.P.F. Evolution Master BUILD UP Brush KIT SIZE 4 build up brush

Item Order No: 2002 M.P.F. Evolution Master BUILD UP Brush KIT SIZE 6 build up brush

Item Order No: 2003 M.P.F. Evolution Master BUILD UP Brush KIT SIZE 8 build up brush



EVOLUTION MASTER BRUSHES



Item Order No: 2004 M.P.F. Evolution Master BUILD UP BRUSH size 4

Item Order No: 2005 M.P.F. Evolution Master BUILD UP BRUSH size 6

Item Order No: 2006 M.P.F. Evolution Master BUILD UP BRUSH size 8

Item Order No: 2010 M.P.F. Evolution Master Stain brush

Item Order No: 2011 M.P.F. Evolution Master Glaze Brush

Item Order No: 2012 M.P.F. Evolution Master Contour Brush



Item Order No: 2007
Replacement brush tip size 4

Item Order No: 2008
Replacement brush tip size 6

Item Order No: 2009
Replacement brush tip size 8

Item Order No: 2013
Replacement Stain brush tip

Item Order No: 2014
Replacement Glaze brush tip

Item Order No: 2015
Replacement Contour brush tip

REPLACEMENT BRUSH TIPS

EVOLUTION MASTER CERAMIC SCULPTING KIT



The Evolution ceramic sculpting kit, this superb range of stainless steel blades have been designed so that the dental ceramist can carve into semi dry ceramic, their fine stainless steel blades are wonderful for carving and shaping semi dry ceramic occlusal areas, these blades can be used to shape your incisal build ups to allow for mamelons and special effects. All blades are replaceable!

We believe this exquisite new line is something that every technician will be proud to work with.

“Be a part of the M.P.F. Evolution; speak to your dealer about this new innovative line of products”.

Item Order No: 2016 M.P.F. Evolution Master Ceramic Sculpting instrument kit



EVOLUTION SCULPTING INSTRUMENTS



Item Order No: 2017 M.P.F. Evolution Master sculpting blade Angle tip



Item Order No: 2018 M.P.F. Evolution Master sculpting blade Javelin tip



Item Order No: 2019 M.P.F. Evolution Master sculpting blade Razor tip



Item Order No: 2020 M.P.F. Evolution Master sculpting blade Flat tip

REPLACEMENT BLADES



Item Order No: 2021
Replacement sculpting blade
Angle tip



Item Order No: 2022
Replacement sculpting blade
Javelin tip



Item Order No: 2023
Replacement sculpting blade
Razor tip



Item Order No: 2024
Replacement sculpting blade
Flat tip

total control

Wax and Die
Management



Wax & Die Management

Wax & Die management is a very important key part of the restoration process in the dental laboratory. We at the M.P.F. Brush Company have taken great care in creating this range of products to allow the technician to have the optimum result from his labor.

Our entire wax and die management range is tailored to achieve precise results when going through the manufacturing process. High quality ash free waxes, die spacers, universal DS3 separator, DH1 die hardener, freeze instant blockout material, debubb for all your debubbler needs, and our exquisite Titanium handled Waxers for comfort and reliability, all comes to you in this exciting new range of products.

Precise hand tools & products designed by technicians for the technician!!!!

Our wax & die products give optimum results
90%
less work for you

PRECISION DIE SPACER



Item Order No: 700



Item Order No: 701



Item Order No: 702



Item Order No: 703

Exhibiting smooth feel and flow, the M.P.F. Precision Die spacer series are a fast drying die spacer and sealer which are characterized by superb thixotropic consistency, uniform homogeneous layer thickness, and precise control of varnish distribution without running or clumping.

Advanced chemistry allows for a strong and compatible bond to virtually any surface, while offering outstanding resistance against pressurized steam, scratching and chipping.

Each coat creates 12 microns of space.

Fast 15-second set time between coats.

M.P.F. Precision Die spacer come in 4 popular colors, gold, blue, silver and grey.



Key Benefits

- Offers accurate die cement spacing.
- Excellent flow properties for precise application.
- Creates a hard, dense surface.
- Fast 15-second drying time.



The ultimate solution for blocking out undercuts on dies.



Item Order No: 706 Base



Item Order No: 707 Accelerator

A full bodied, controlled viscosity material used for blocking out undercuts, filling air bubbles, and building cusp tips on dies prior to wax up.

Offering excellent contouring and shaping characteristics, Freeze is applied to the area to be blocked out or based up. The applied base is then sprayed with Freeze Accelerator for instant curing. The result: A hard, smooth, glass-like surface true to the surface contour of the die. Freeze Accelerator can also be used to instant cure all cyanoacrylate adhesives on the market.

Key Benefits

- Simplify die repair procedures.
- Join wax, metal or porcelain.
- Position attachment and implants.
- Burns out clean.

Caution: Freeze base material should not be positioned in the same place on the bench, as the accelerator fumes can inhibit the base material and cause it to harden.

FREEZE BLOCKOUT




1. Spray a little accelerator into the region you want to fill with base material.
2. Then apply the base material carefully.
3. Spray accelerator directly onto the base material and allow it to set, for 20sec.
4. After the base material is hard, it can also be cut back with a fine grade rotary instrument for better contouring.

Ti WAXER




Ti Waxer Blue Code handle
Item Order No: TW-02

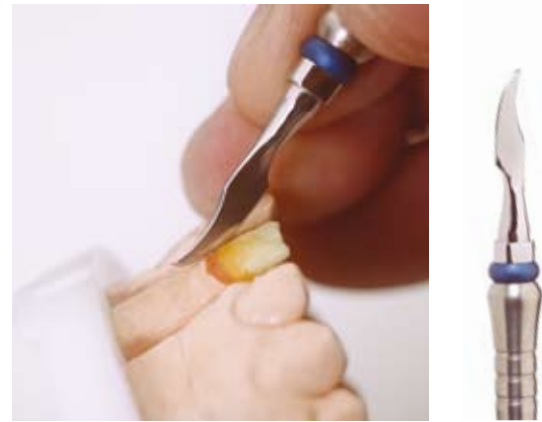

Ti Waxer Red Code handle
Item Order No: TW-03



Ti Waxer Purple Code handle
Item Order No: TW-04

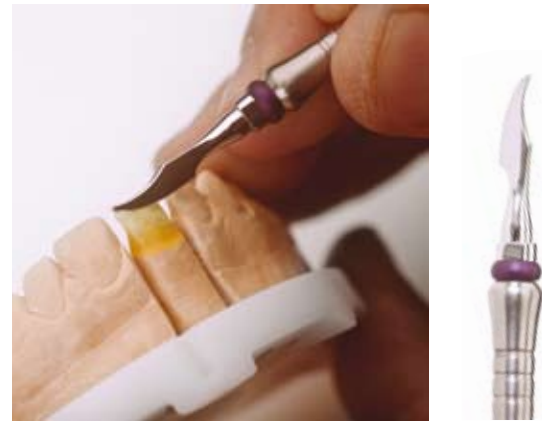

The M.P.F. Titanium handled waxing instruments are a latest innovation from the M.P.F. Brush Company. These elegant wax carvers have been designed with comfort in mind, the long hours we spend waxing up delicate crowns can be tiresome, our light weight handle and ergonomic design make the technicians life a little easier.

With our double-sided Ti waxer the technician has instant access to two different carving instruments, with a quick turnaround of the handle the technician has the ability to use the instrument for several applications, carving, connecting!
The stainless steel blades are replaceable and each handle is color coded so the technician can identify which handle he needs at a quick glance.



The Ti Waxers are sold in a 3 in 1 kit box, or they can be bought individually.

M.P.F Ti Waxer 3 in 1 kit Item Order No: TW-01

The M.P.F. Ti Waxer can be used with a regular Bunsen flame and also the latest "No Flame" systems.



All components of the T Waxer can be purchased individually.

M.P.F. Ti Waxer, Blue Code, replacement large Lecron blade Item Order No: TW-05

M.P.F. Ti Waxer, Blue Code, replacement tapered Spade blade Item Order No: TW-06

M.P.F. Ti Waxer, Red Code, replacement Spoon Item Order No: TW-07

M.P.F. Ti Waxer, Red Code, replacement large Dropper Item Order No: TW-08

M.P.F. Ti Waxer, Purple Code, replacement small Lecron Item Order No: TW-09

Ti Waxer, Purple Code, replacement small Dropper Item Order No: TW-10



WAX FORM MODELING & SCULPTING - WAX



Item Order No: 714
Sculpting Wax | Beige



Item Order No: 715
Sculpting Wax | Blue



Item Order No: 716
Modeling Wax | Green



Item Order No: 717
Modeling Wax | Ivory



For creative, inspired sculpturing of dental morphology.

Technical Application

Primarily formulated for occlusal anatomy, and sculpting full crowns, inlays and onlays with ease. Suitable for pressable ceramics.

Working Properties

A pure, ash-free European wax composition with excellent flow properties and low contraction for modeling the exterior tooth shape. Excellent carvability without compromising strength. Exhibits control and accuracy.



Key Benefits

- Offers minimal shrinkage.
- Excellent sharp carving characteristics.
- Ash-Free wax.
- Suitable for pressable ceramics.

DEBUBB DEBUBBLIZER



Item Order No: 708

M.P.F Debubb is an advanced debubblizer clinically designed to reduce surface tension, thus eliminating bubble formation when pouring dental model casts from impressions and investing wax patterns. Innovative trigger mist™ delivery system offers ergonomic comfort and sprays an ultra-fine mist to ensure even coverage. Safe for use on all types of impression materials, including alginates, wax and hydrocolloids. M.P.F Debubb is non-contaminating and alcohol free.



Key Benefits

- Reduces surface tension on wax and impressions.
- Allows gypsum and investment to flow smoothly into all areas.
- Ensures bubble free models and smooth precision castings.
- Guarantees consistent, accurate results.

PEGFIX REFRACTORY MATERIAL



Item Order No: 704

Peg Fix is a second generation, advanced high-heat resistant refractory material used for custom ceramic firing stabilizing. The material can be easily shaped and used with any type of porcelain.

Available for a variety of uses in the laboratory, the material is stable up to 1300°C (2375°F), eliminates firing distortion and contamination. Pegfix offers a freeze stable solution for hassle free shipping and extended shelf life for added value. Packaged in 12cc Syringes.



Key Benefits

- Fill inside of crowns for extra stabilization during ceramic firing cycles.
- Repair cracked pressable ceramic cores or all-ceramic margins.
- Use for stabilization during soldering procedures.
- Use to hold platinum foil in place during cast repairs.

DS3 DIE SEPARATOR



Item Order No: 710 DS3 Pen Applicator



Item Order No: 712 DS3 Bottle & Brush

One separator, three applications. The new DS3 Die Separator is an advanced universal separating medium used for the clean separation of Porcelain Margins, Composite Restorations and Wax Copings from the die. The thin liquid evaporates quickly, leaving a non-oily surface with no layer build-up. The material is available in a felt pen applicator, 15ml brush bottle.



Key Benefits

- 1 advanced separator, 3 technical applications.
- Non-oily with no layer build-up.
- Non-contaminating.
- Fast 10-second drying time.



DH1 DIE HARDENER & SEALER



Item Order No: 709 DH1 Pen Applicator



Item Order No: 711 DH1 Bottle & Brush

Seals and strengthens all die materials enabling accurate die preparation.

A paint-on acrylic resin for conditioning, sealing and hardening gypsum dies and models. Water proofs and seals stone surface without layer build-up on the die. Advanced material is absorbed into the stone and by capillary action, seals the stone surface while providing zero micron buildup. Protects margins and reduces abrasion. DH1 is available in a felt-pen applicator, as well as a 15ml bottle brush.


Key Benefits

- Strengthens dies and hardens models.
- Protects margins from chipping and abrasion.
- Zero microns ensure zero layer buildup.
- Improves margin fit.



concentric
& reliable

Rotary
Instruments

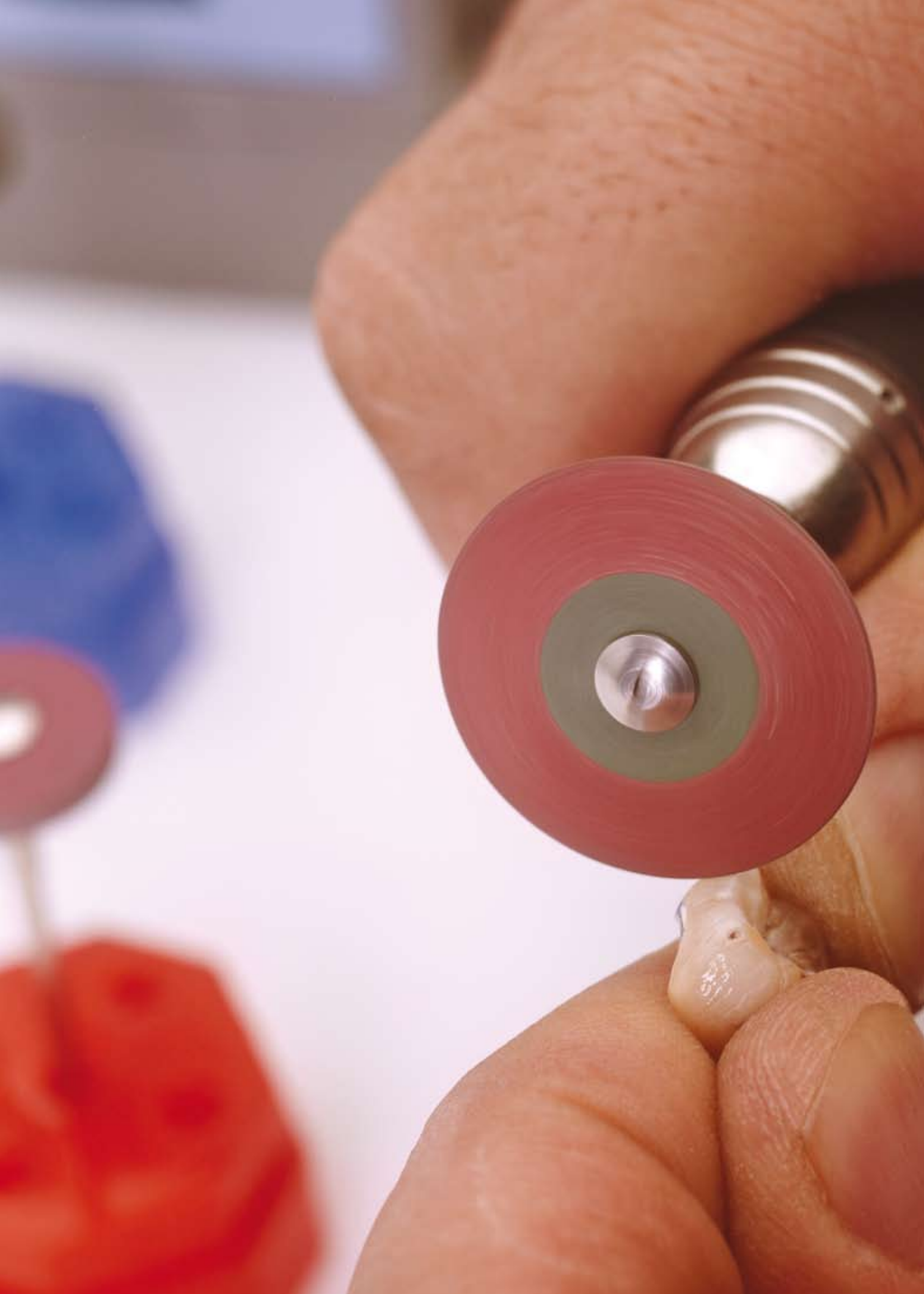

Rotary Instruments

The Dental Ceramic Industry is changing constantly.

As we have seen a steady trend towards the Cad Cam Industry utilizing such wonderful materials such as Alumina and Zirconia, this trend has led us to focus our Rotary Instrument production in this direction.

We believe that this trend will continue to expand and move forward year after year, so the production of high quality abrasive tools for these new materials is our future.

Our range of Diamond Impregnated Rubber Discs, Zirco-Cera and our GSS Green Stone Series are constantly finding their foundations within this competitive Dental Abrasives market around the world today.

Due to zirconia's hardness and sensitivity to heat, there was a great need for a rotary instrument that could cut and abrade this material with a minimum of heat generation as possible. **M.P.F.** has created it!! Our latest Innovation the **Z maX** Diamond Impregnated wheels.

We have created a rotary instrument with a very high content of fine diamond particle which has been bonded with special binders that hold all the elements together. This creates a superb Zirconia cutting instrument with exceptional life span, **2-3 times more life than any product on the market today.** With our high content of diamond particle our product can cut and abrade with very little pressure needed from the hand piece, creating less heat generation to the zirconia substructure.

We at the M.P.F. Brush Company believe that a high quality product together with a fair price is the way forward to a successful future.

Try our products and let us know if we have achieved this goal!

Increased diamond content means

2 TO 3
times longer lifespan
than any other product
on the market today

RAINBOW CUTTERS

These diamond impregnated rubber wheels are exceptional quality, they are to be used for all types of dental ceramic polishing and can be used to abrade and polish zirconia.

They come in four grades Coarse Red, Medium Blue, Grey Fine and Yellow Super Fine for a luster shine.

The two tone green color coding system can allow you to distinguish from other products sold in the market.

We advise that you stick to the RPM speeds to get the optimum result and by doing so your product will have an extended life span.



Red Coarse **Item Order No: RC-0001** RPM 13000 - 15000



Red coarse rubber disc a highly impregnated rubber disc for minimal reduction on zirconia marginal regions and all types of ceramic materials with great ease, this range of rubber wheels has exceptional long life properties with all ceramic materials.
RPM 13000-15000




Blue Medium **Item Order No: RC-0002** RPM 13000 - 15000



Blue medium rubber disc a highly impregnated rubber disc for minimal reduction on zirconia marginal regions and for contouring all types of ceramic materials with great ease, this range of rubber wheels has exceptional long life properties with all ceramic materials.
RPM 13000-15000




Grey Fine **Item Order No: RC-0003** RPM 10000 - 13000

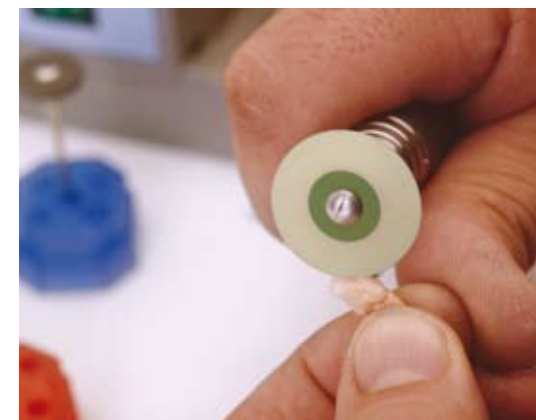

Grey Fine rubber disc a highly impregnated rubber disc for the removal of minor scratches on all types of ceramic materials with great ease, this range of rubber wheels has exceptional long life properties with all ceramic materials.
RPM 10000-13000



Yellow Super Fine **Item Order No: RC-0004** RPM 8000 - 10000

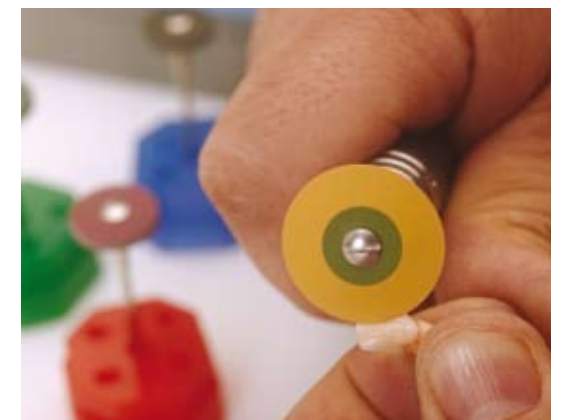


Yellow Superfine rubber disc a highly impregnated rubber disc for the final high shine procedure on all types of ceramic and zirconia materials, this range of rubber wheels has exceptional long life properties with all ceramic materials.
RPM 8000-10000

These 4 grades of premium diamond abrasive wheels ensure a perfectly sealed ceramic or zirconia restoration. By using our M.P.F. Diamond Paste you can be sure to obtain a highly polished, finished restoration

Z max THE ULTIMATE ZIRCONIA ABRASIVE WHEEL

Z Max, highly impregnated diamond wheels and discs, to abrade Zirconia.

The Z max series comes in 9 different shapes:

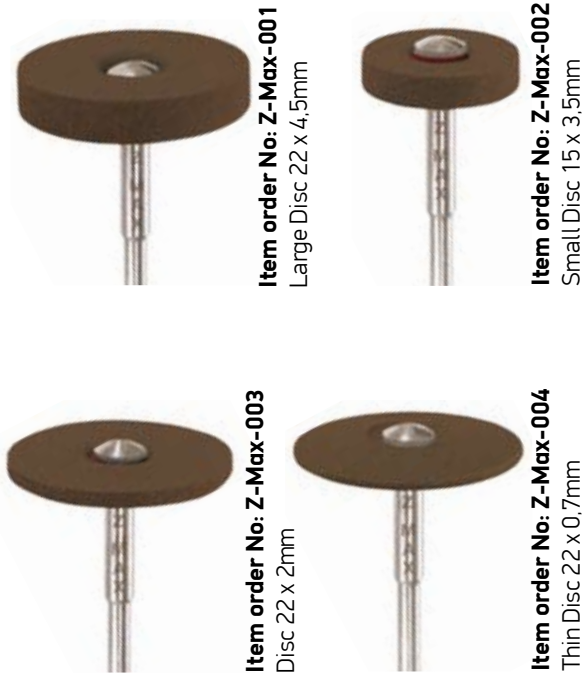

Zirconia has become very popular in the dental industry in the past few years, and its growing day by day.

Due to zirconia's hardness and sensitivity to heat, there was a great need for a rotary instrument that could cut and abrade this material with a minimum of heat generation as possible.

M.P.F. has created it!!

We have created a rotary instrument with a very high content of fine diamond particle which has been bonded with special binders that hold all the elements together, creating a superb Zirconia cutting instrument, with exceptional life span, 2-3 times more life than any product on the market today.


With our high content of diamond particle our product can cut and abrade with very little pressure from the hand piece, thus much less heat generation to the zirconia substructure. Utilizing stainless iron mandrels for extra strength we guarantee a very high end looking and functioning product.

We have laser printed onto our mandrels for specific marking & branding for this exclusive product line.

DIAMOND IMPREGNATED SILICON RUBBER WHEELS

Zirco-Cera Diamond Polisher

Used to cut and polish Zirconia and all types of ceramic material.

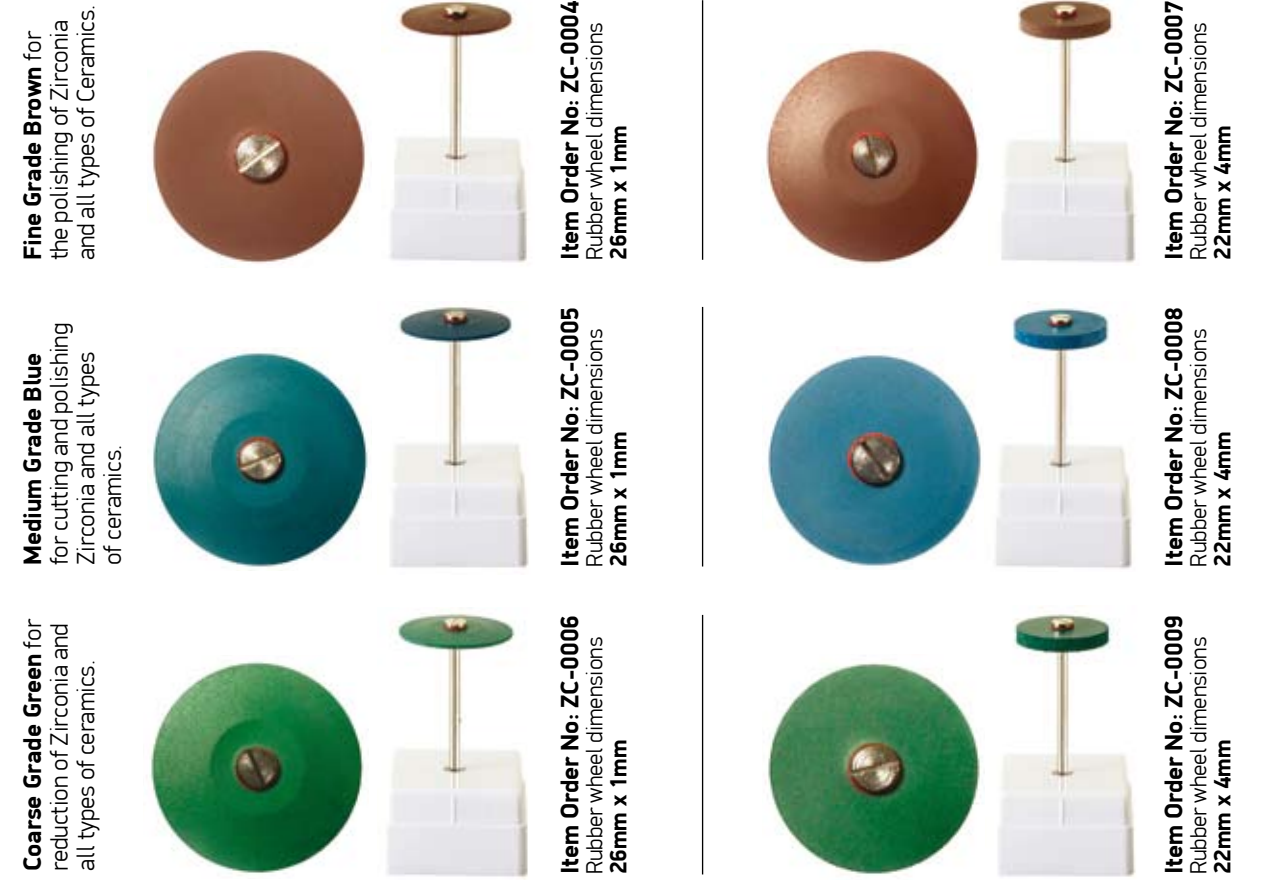

Specifically for the cutting and polishing of all dental ceramic materials.



RUBBER WHEELS & POINTS

Granular Grade: COARSE

Packaging (piece/box):
F & H Series = 100
K Series = 50, 12, 6
Piece / Box

ID: Products ID
Diameter / Height
Max R.P.M.

Unmounted F Series



FB006221	FC008221	FA001061
22/3	22/4	6.5/24
10000	10000	10000

Unmounted H Series



HA001061	HB00221	HC00221
6/22	22/3	22/4
20000	10000	10000

Granular Grade: MEDIUM

Packaging (piece/box):
F & H Series = 100
K Series = 50, 12, 6
Piece / Box

ID: Products ID
Diameter / Height
Max R.P.M.



FB006222	FC008222	FA001062
22/3	22/4	6.5/24
10000	10000	10000



HA001062	HB00222	HC00222
6/22	22/3	22/4
20000	10000	10000

Granular Grade: FINE

Packaging (piece/box):
F & H Series = 100
K Series = 50, 12, 6
Piece / Box

ID: Products ID
Diameter / Height
Max R.P.M.



FB006223	FC008223	FA001063
22/3	22/4	6.5/24
10000	10000	10000



HA001063	HB00223	HC00223
6/22	22/3	22/4
20000	10000	10000

Silicon rubber wheels and points designed to cut polish, semi precious and precious alloys. Specifically for dental laboratory use.

Silicon rubber wheels and points designed to cut and polish dental ceramics. Specifically for dental laboratory use.



Mounted K Series



KA222131	KA226101	KA228151
13/20	10/24	15/18
8000	10000	8000



KA222132	KA226102	KA228152
13/20	10/24	15/18
8000	10000	8000



KA222133	KA226103	KA228153
13/20	10/24	15/18
8000	10000	8000

Designed specifically to cut and polish acrylic resin base plates, denture base. Specifically for dental laboratory use.



Mounted Mandrels



The **M.P.F. Premium Stainless Steel Mandrels** are made from a high quality SS Grade Steel.

Shank size 2.35mm



DIAMOND IMPREGNATED GREEN STONE SERIES

Zirco-Cera Abrasive Wheels

"GSSeries have the ability to cut all types of Zirconia and Ceramic material with great ease... The surface is left with a smooth, clean finish."

The GSSeries comes in 9 different shapes:



Item order No: GSS-005
Small Barrel 4,8 x 13mm



Item order No: GSS-006
Large Barrel 6,5 x 13mm



Item order No: GSS-007
Cone 3,5 x 11mm



Item order No: GSS-008
Inverted Cone Large 12 x 6mm



Item order No: GSS-009
Inverted Cone Small 6 x 8mm



Item order No: GSS-004
Thin Disc 22 x 0,7mm



Item order No: GSS-003
Disc 22 x 2mm



Item order No: GSS-002
Small Disc 15 x 3,5mm



Item order No: GSS-001
Large Disc 22 x 4,5mm

The GSSeries is a diamond impregnated rotary instrument of the highest quality, boasting a very high content of industrial diamond, bound together with a special ceramic binder. The GSSeries rotary instruments have extraordinary abrasive ability, they can cut through Zirconia with a minimal amount of pressure, this alleviates the chances of micro fractures within the Zirconia substructure.

GSSeries have the ability to cut all types of ceramic material with great ease and due to the fine grain. Industrial diamond and special binder, the ceramic surface is left with a smooth, clean finish. Due to the composition of the diamond to binder ratio, the GSS Series has a very high service life, ergonomically it is by far the best solution for the ceramist today!



Item order No: 037
Disc dimensions: 25 x 0,6 x 1,8mm



Item order No: 038
Disc dimensions: 38 x 0,6 x 1,8mm



Item order No: 039
Disc dimensions: 22 x 0,2 x 1,8mm



ALUMINIUM OXIDE ABRASIVE DISC

The durability of our discs second to none!

This range of Aluminum Oxide Abrasive Discs are of the highest quality in the market today. Our discs have the ability to work on all types of alloy, from Gold, to the hardest non precious alloys the market has to offer.

No metal is too hard for the M.P.F. Aluminum Oxide Abrasive Disc, a must for the partial denture laboratory!

The 25mm disc is a work horse, it has the ability to cut units from the sprue, but also trim and prepare your metallic surface with very little pressure.

38mm disc are great for cutting units away from their sprue!

Due to the binder that we use to make our discs with, we have the ability to use our discs in ways that are not conventional to the normal disc, left to right movements are no problem for the M.P.F. Abrasive Disc, the binder allows us to have this flexibility, without breakage.

Our super fine 0.2 - 0.25mm disc gives us several options:

1. We have the ability to cut and separate all types of dental ceramic, without it chipping and cracking the ceramic surface.
2. We can also cut off units from our High Noble content alloys with great ease, and with very little waste, due to the thin 0.2 disc size.

POLISHING BRUSHES FOR CERAMIC & METAL



Item order No: 040
Sold in quantities of 100, 50 or 25 pieces per pack.

Our polishing brushes are made from natural hair and polish all metallic surfaces extremely well. If they are used at the correct RPM they can have a long service life.



Item order No: 041
Sold in quantities of 100, 50 or 25 pieces per pack.

Our polishing mop is the perfect hardness and thickness for polishing ceramic surface when combined with the M.P.F. Diamond polishing paste, you are certain to get a brilliant high shine surface.

TUNGSTEN CARBIDE BUR

F.G. DIAMOND BUR

Granular Grade: COARSE

Packaging (piece/box): 1



ID: Products ID
Head Diameter
Color code: Green

Granular Grade: SUPER COARSE

Packaging (piece/box): 1



ID: Products ID
Head Diameter
Color code: Black

Granular Grade: FINE

Packaging (piece/box): 1



ID: Products ID
Head Diameter
Color code: Red

Granular Grade: SUPER FINE

Packaging (piece/box): 1



ID: Products ID
Head Diameter
Color code: Yellow



Medium to Coarse Grade: Blue



ISO No. 170/018 199/016 197/016 197/018 160/014 160/012 200/012 289/014 140/012 141/014 110/014
ID: Products ID TF-23 TR-12 TR-21 TR-24 TC-21 TC-26 SO-20 SO-21 SR-11 SR-12 SF-21

Medium to Coarse Grade: Blue



ISO No. 257/018 111/012 111/014 234/014 237/010 237/018 010/012 010/018 001/012 001/016 068/042
ID: Products ID FO-32 SF-11 SF-12 EX-31 EX-41 EX-20 SI-46 SI-49 BR-46 BR-40 WR-13

Medium to Fine Grade: Red



ISO No. 173/016 199/016 199/018 160/016 298/014 298/016 237/021 196/020 197/012 254/017 108/008 138/007
ID: Products ID TF-12F TR-11F TR-26F TC-11F FO-21F FO-22F EX-21F CR-11F CR-21F FO-30F CD-58F CD-52F

Our specially selected diamond friction grip burs come in 22 different shapes in our medium to coarse grade, that are color coded in blue and 12 different shapes in our medium to fine grain which are color coded in red.

Each color and shape is packaged in plastic boxes of 5 pieces per box.

purest form

Techlight



LED Technology

The Future of Lighting

The M.P.F. Tech Light boasts the latest lighting technology that is making rapid changes in the world today. This unique lighting technology offers lighting in its purest form. We at the M.P.F. Brush Company have designed this light in such a way, that it is not only light source, it's much more!

- + 5,500 Kelvin of pure daylight.
- + Telescopic pole for height adjustment.
- + Multifunctional arms to adjust the light into the position you require on your lab bench.
- + Power outlet available to plug in a wax pot, micro motor, under 500 watt.
- + 5 Super Power LED chip lamps to light up your bench. Separate switch for the 5 lamps.
- + 1 Super Power LED chip lamp. This pilot light is designed to give off a very strong beam of light and has 180 Degrees maneuverability. It can be used from the front or from the back of the light housing. This pilot light has been designed to allow you spectacular visual of the deepest crown or any problem area that needs to be investigated, this light makes a huge difference in tight spots that are very difficult to see. Separate switch for the pilot light.
- + The light housing has the ability to rotate also at 180 Degrees, giving this light the ability to make in house shade taking, much easier.
- + The light structure is made from a combination of Stainless Steel and Aluminum that has been color anodized to protect the metal surface.

Our revolutionary
techlight provides
5500kl
Pure daylight
Perfect clarity!

..... We like to thank **Eladio Gonzalez Gomez**
of "**Fusion Taller Dental**" Mexico
E tallerfusion@prodigy.net.mx

TECHLIGHT



Pilot Light Manouverability

The pilot light is a unique design, which is a must for dental laboratory work. This super LED lamp shoots out a very strong beam of light, this is essential for working in very tight dark areas such as the depths of very lengthy crowns, intricate attachment work, complicated implant screws and abutments. This light makes these problem areas look like "daylight". The pilot light also has 180 degree manouverability. It can be used either from the back of the housing or from the front.

Light Switch

The light has two separate switches. The top switch is for the pilot light and the bottom is for the 5 lamps located in the housing.

Light Housing Ability

The complete housing can turn 180 degrees in both directions. This manouverability makes in-house shade taking much easier.

N.B. Patient should wear protective glasses, as the super LED lamps are very strong.



M.P.F. Techlight
Replacement Bulbs
Item Order No: 057

1. Pilot light Super LED chip light.
2. 5 Super LED chip lights in a row.
3. Two separate switches for the pilot light and the 5 lights in the housing.
4. Multifunctional arms allow you to maneuver the light into the position you require.
5. Telescopic column can be adjusted to the height required.
6. This brace is used to tighten the telescopic arm into position, by simply tightening in clockwise direction, the column will stabilize.
7. AC power outlet can support any appliance up to 500 watt (e.g. Wax pot, micromotor).


expertise
& excellence

Pilot Lab



Nondas
Vlachopoulos

From the very beginning of our endeavour with the M.P.F. spring activated ceramic brush, we have had the input and support of one of the finest dental technicians we have come across in 30 years in the industry. His name is **Nondas Vlachopoulos**. Nondas is a Canadian born Greek, who has studied in Europe and has excelled in his field in a big way. He has held seminars in Europe and the USA. His expertise and direct input has been very important to the development of our M.P.F. product range. We would like to take this opportunity to thank Nondas for all his support and thank him for the wonderful photos of the crowns that you see displayed throughout our catalogue, **simply amazing!**

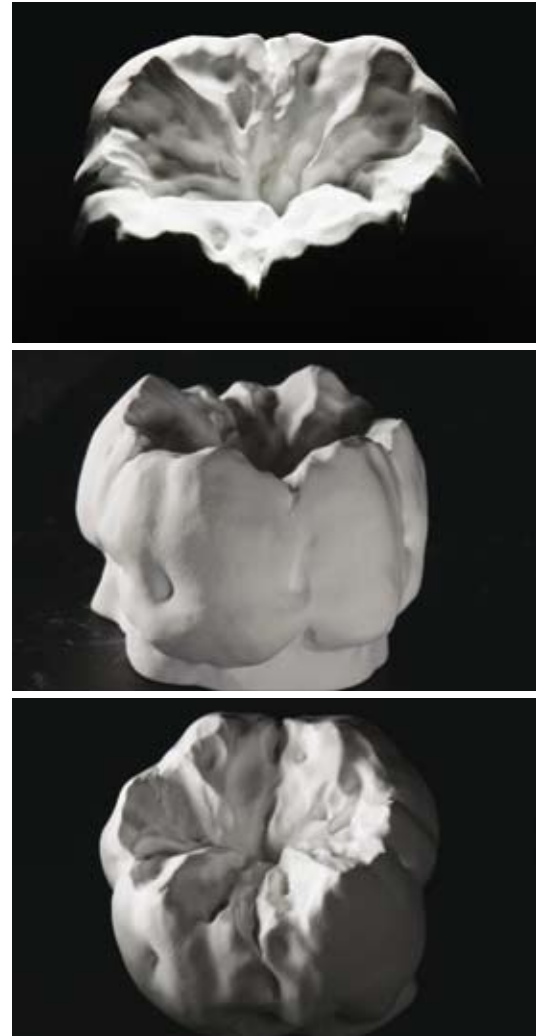

If you would like to have more information on this exceptional dental technician, you can contact him at: **Aesthetic Lab**

Dental Technician Nondas Vlachopoulos
Pleiadon Str. 45, Nea Kifisia, 145 64, Athens, Greece
T +30 210 6250287 | +30 210 6250288
E nondasdt@hol.gr **Facebook:** Aestheticlab
www.aestheticlab.gr

We ask the best
what they think
giving us over
30
years of expertise
at our finger tips



Peter
Pizzi

Over the past 18 months we at the M.P.F. Brush Company have had the pleasure of working very closely with top dental clinician, **Peter Pizzi**. As a motivational dental educator, Peter has found an easy transition into the lecture/clinician circuit. His personal appreciation and unique lecture style has helped to expand the lecture concept and create an informative, comfortable environment. Peter's expertise on all phases of crown and bridge, porcelain, implantology, muscle function, mandibular physiology and esthetic communication have made him a source of knowledge and motivation for his peers. These past 18 months have been a very exciting time for both Peter and the M.P.F. Brush Company. From this wonderful collaboration we've seen the birth of a very high end line of dental ceramist hand tools under the name **M.P.F. Evolution**. With Peter's hands on approach to pretty much everything he is involved in, he has brought to life this unique range of precision instruments. Peter lectures to dentists and technicians throughout the United States, Europe, South America and Asia on several different topics. **"We at the M.P.F. Brush Company are very fortunate and proud to have Peter Pizzi as our pilot lab and technical guidance counselor in the USA."**



For information on this very talented dental clinician you can visit his web site at **www.pizzidentalstudio.com** or contact Peter at: **Pizzi Dental Studio**

4038 Victory Blvd, Staten Island, NY 10314
T 718-477-3771 E pizzi7@aol.com



M.P.F.
BRUSH COMPANY

44 Alekou Panagouli St. Rhodes 85100 Greece
T | F +30 22410 37336 **E** info@mpfbrush.com
Facebook www.facebook.com/mpfbrushcompany
Skype cosdent **www.mpfbrush.com**

